

Pro-Bario Active

Bariatric Community Bed



The Pro-Bario Active bed is designed for plus-size users living at home with complex needs.


With a safe working load of 380kg (60st) and a wide sleep deck, it provides a comfortable and protective platform for the user.

Auto-contour and auto-regression help to prevent the user sliding down the bed as they sit up, reducing skin shear and friction.


A range of compatible bariatric mattresses and accessories including chairs, cushions, walking aids, hoists and a commode/shower chair is available from Medstrom Select.



A safe and practical home care solution for plus-size users and their carers

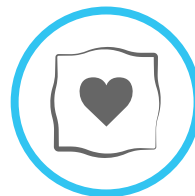

Wide platform provides user comfort and space for caregivers to work safely

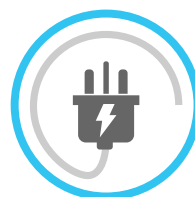

Electric adjustment reduces manual handling risks

Additional Benefits

The Pro-Bario Active community bed also provides these additional benefits as standard:

- **Auto-contour and auto-regression reduce sliding down the bed, reducing shear and friction.**
- **Top height allows tall caregivers to work safely when providing personal care.**
- **Integral side rails, fully compliant to IEC 60601-2-52, for user safety.**
- **Powder-coated steel frame is durable and easy to clean.**
- **Consists of four parts, which are carried into the home separately and safely assembled inside.**
- **Range of compatible mattresses and accessories available.**



Technical Specification

Safe Working Load (SWL)	380kg / 60st
Max Patient Weight (MPW)	350kg / 55st
Height Range	30cm - 70cm
Sleep Deck Dimensions	120cm x 200cm
Bed Frame Dimensions	136cm x 217cm
Backrest Angle:	to 70°
Thigh Section Angle:	to 30°
Calf Section Angle:	to 20°
Trendelenburg / Reverse Trendelenburg	10° to -10°